

# JumpStart

ANTIMICROBIAL WOUND DRESSING



**Inspired** by the body.  
**Powered** by electricity.  
**Energized** by results.®

**Arthrex**®  
**Vet Systems**





# JumpStart™


## Antimicrobial Wound Dressing

Used in multiple applications, JumpStart provides sustained broad-spectrum antimicrobial efficacy and an optimal environment for wound healing.


# Inspired by the body.

## Electricity Is Essential to Healing


The body naturally creates and uses electrical energy to promote healing.<sup>1</sup>

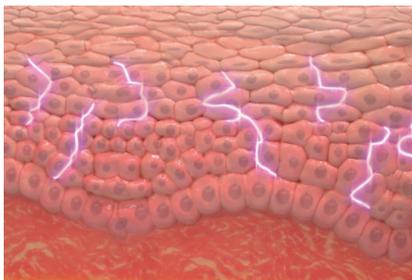

An electric potential exists across intact skin.<sup>2-4</sup>

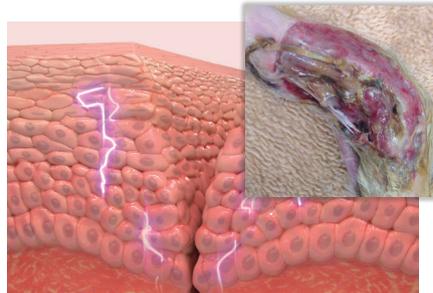

When skin is wounded, a change in electric potential occurs. This stimulus is the earliest guidance signal to initiate cell migration and re-epithelialization, and is essential to wound healing.<sup>5</sup>

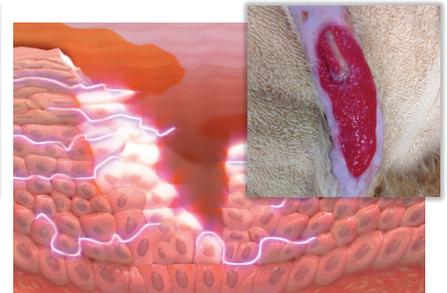


Designed to mimic the body's natural physiologic electric currents, JumpStart's embedded microcell batteries generate microcurrents in the presence of a conductive medium<sup>6</sup> to harness the power of electricity and support the body's natural electrical healing process.<sup>7</sup>

## Versatile Application

JumpStart is used for a broad variety of superficial and full-thickness wounds:

- Surgical incisions
- Donor and/or recipient graft sites
- First- and second-degree burns
- Non-healing wounds
- Traumatic wounds
- Cuts and abrasions
- Punctures and lacerations
- Irritations and abscesses
- Proud flesh

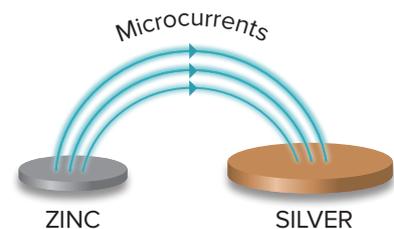

JumpStart antimicrobial dressing is powered by patented V.Dox™ technology – employs moisture-activated microcell batteries that wirelessly generate microcurrents.

# JumpStart™ Dressings Product Features


## Advancing the Science of Animal Wound Care

JumpStart antimicrobial wound dressings use electricity to provide broad-spectrum antimicrobial protection for incisions and wounds in both small and large animals.

- Conforms easily to body contours for patient comfort
- Multiple sizes available to fit surgical and wound sites
- Can be cut to fit under secondary dressings
- Designed for multiday use for long-lasting protection



Proprietary islands of elemental silver and zinc form microcell batteries that generate electricity in the presence of moisture

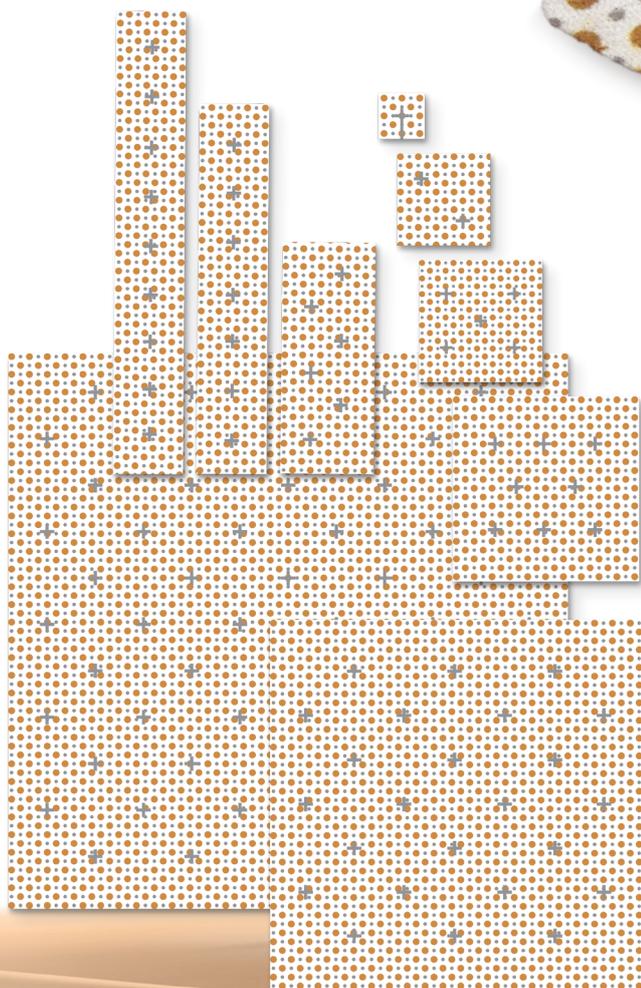


# Powered by electricity.

## Augmenting the Body's Natural Healing Process

---

JumpStart™ Wound Dressings are an entirely new generation of incision and wound care solutions.

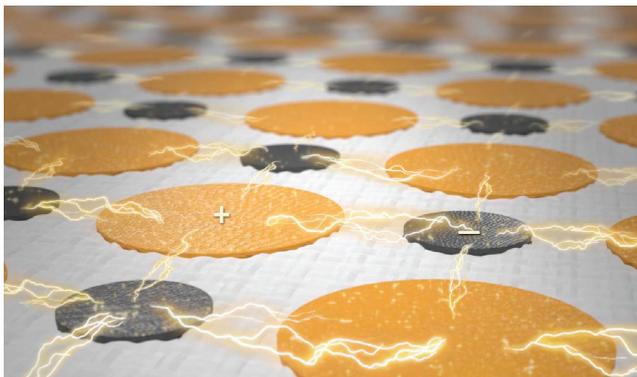

### Advanced Microcurrent Technology

---

- Generates physiologic levels of microcurrent (2-10  $\mu$ A) that are known to be necessary for healing<sup>1</sup>
- Flexible and portable with no need for external power source



### Comfortable and Conformable

---

- Conforms to wound for maximum mobility and comfort
- Easily cut to fit
- Available in a variety of configurations, including JumpStart composite dressing and JumpStart non-adhesive single-layer dressings



### Broad-Spectrum Antimicrobial Efficacy

---


JumpStart's antimicrobial protection<sup>8-10</sup> creates an optimal environment for wound healing.<sup>7</sup>

- Provides effective and sustained antimicrobial efficacy<sup>11</sup>
- Kills antibiotic-resistant and biofilm-forming pathogens<sup>8-10</sup>
- No silver release into the bloodstream<sup>6</sup>


# Energized by results.

## Optimizing Outcomes

Studies have shown JumpStart™ reduces the risk of infection<sup>9,10</sup> and promotes the healing process<sup>7</sup> to optimize outcomes.

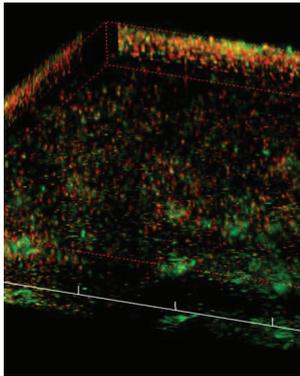

JumpStart dressing

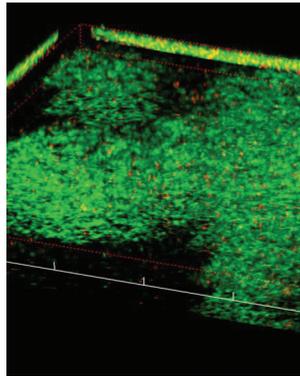

Silver only

## Reduce Risk of Infection



- Effective, sustained broad-spectrum antimicrobial efficacy for up to 7 days<sup>11</sup>
- Prevents bacterial growth<sup>10</sup>
- Electricidal antimicrobial efficacy vs silver alone<sup>10</sup>

Live/dead fluorescence staining demonstrates bacterial killing of *P. Aeruginosa* within JumpStart dressing compared to a standard silver-based dressing at 24 hours. Red = dead, Green = alive



## Optimal Environment for Wound Healing

- Keratinocyte migration vs silver observed at 9 hours, in vitro scratch assay<sup>7</sup>

JumpStart Dressing

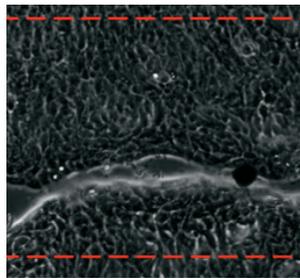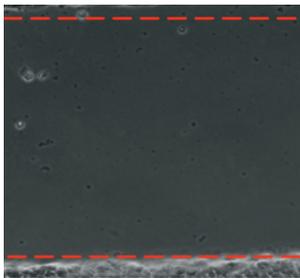

Silver Only

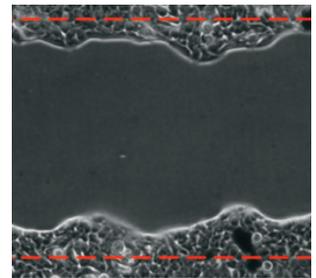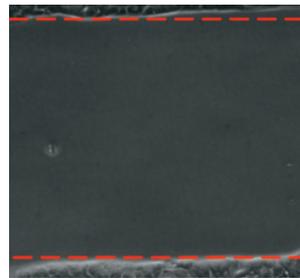

Initial wound



One month later

## Promote Healing<sup>12-15</sup>

- Improved scar appearance vs standard dressings<sup>13,15</sup>
- 45% shortened wound healing time in a clinical study<sup>14</sup>
- Cleared by the US FDA for partial- and full-thickness wounds in humans

# Electrical Antimicrobial Efficacy

JumpStart™ kills a broad spectrum of harmful pathogens, including multidrug-resistant<sup>8</sup> and biofilm-forming bacteria<sup>9,10</sup> to help reduce the risk of infection.

Biofilm-Producing Pathogens	Anti-Biofilm Efficacy
<i>Acinetobacter baumannii</i>	4
<i>Corynebacterium amycolatum</i>	4
<i>Enterobacter aerogenes</i>	4
<i>Enterococcus faecalis</i>	4
<i>Escherichia coli</i>	4
<i>Klebsiella pneumoniae</i>	4
<i>Pseudomonas aeruginosa</i>	4
<i>Staphylococcus aureus</i>	4
<i>Staphylococcus epidermidis</i>	4
<i>Serratia marcescens</i>	4
Antibiotic-Resistant Pathogens	Bactericidal Efficacy
<i>Klebsiella pneumoniae (ESBL)</i>	4
<i>Pseudomonas aeruginosa (MDR)</i>	4
<i>Staphylococcus aureus (MRSA)</i>	4
Vancomycin-resistant <i>Enterococcus raffinosus</i> (VRE 510)	4
<i>Staphylococcus aureus</i> (NRS1, NRS12, NRS73, NRS116)	4
Vancomycin-resistant <i>Staphylococcus aureus</i> (VRS1, VRS9, VRS11b)	4

Pathogens	Antimicrobial Efficacy
<i>Acinetobacter baumannii</i>	4
<i>Acinetobacter calcoaceticus</i>	4
<i>Aspergillus niger</i>	4
<i>Candida albicans</i>	4
<i>Corynebacterium amycolatum</i>	4
<i>Corynebacterium xerosis</i>	4
<i>Enterobacter aerogenes</i>	4
<i>Enterobacter cloacae</i>	4
<i>Enterococcus faecalis</i>	4
<i>Escherichia coli</i>	4
<i>Klebsiella pneumoniae</i>	4
<i>Propionibacterium acnes</i>	4
<i>Pseudomonas aeruginosa</i>	4
<i>Serratia marcescens</i>	4
<i>Staphylococcus aureus</i>	4
<i>Staphylococcus epidermidis</i>	4
<i>Staphylococcus simulans</i>	4
<i>Streptococcus pneumoniae</i>	4
<i>Trichophyton rubrum</i>	4

The tables show the efficacy of JumpStart dressings against various bacterial and microbial species in vitro, using AATCC Test Method 100-1993.<sup>8</sup>

Some of the strains used for efficacy testing were clinical wound pathogens. In some of the studies more than one strain for each organism was tested. Anti-biofilm efficacy testing was conducted using poloxamer hydrogel and/or colony-drip flow reactor (CDFR) biofilm models.

# Improved Clinical Outcomes



Equine - Initial



Equine - One month



Feline - Initial



Feline - Day 17



Canine - Initial



Canine - Week 1




Canine - Initial



Canine - Week 1

# Contact Layer Instructions for Use

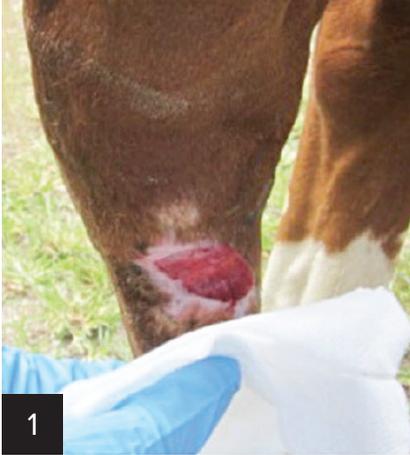

Cleanse the wound area with an appropriate water-based wound cleanser.



If needed, cut to a shape that will extend 1 to 2 inches beyond the wound edge.



Moisten the JumpStart™ dressing with a sterile agent such as saline solution, water, or water-based wound gel (hydrogel).



Apply JumpStart dressing to the wound site, with the dotted side against and directly contacting the wound surface.




Cover JumpStart with two layers of moistened sterile gauze to maintain appropriate moisture levels. Secure in place with elastic adhesive tape.




To ensure even pressure distribution, apply cotton batting to the wounded area and surrounding tissue.



7  
Brown gauze may be used to secure the cotton batting layer.

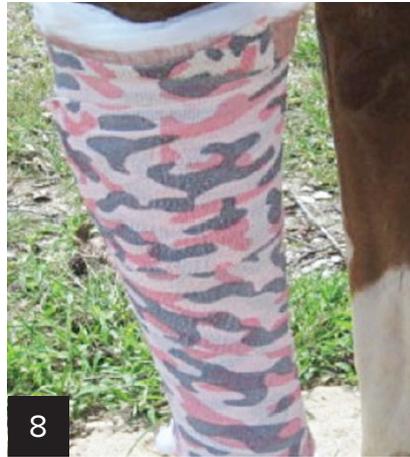

8  
Apply a self-adherent bandage to secure the dressing in place and to extend its longevity.

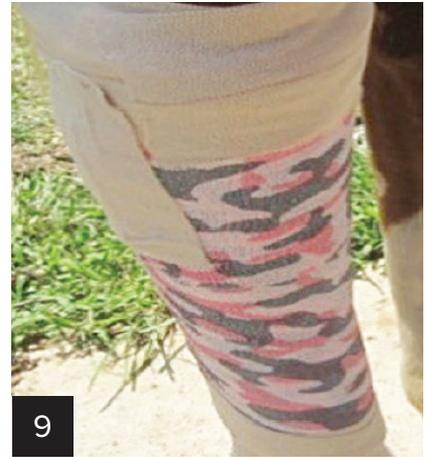

9  
Apply cotton elastic adhesive tape to the dressing edges to secure in place.

For full application instructions and safety information, please refer to JumpStart™ Wound Care with patented V.Dox™ technology Instructions for Veterinary Use.

## References

1. Nuccitelli R. A role for endogenous electric fields in wound healing. *Curr Top Dev Biol.* 2003;58:1-26.
2. McCaig CD, Rajnicek AM, Song B, Zhao M. Controlling cell behavior electrically: current views and future potential. *Physiol Rev.* 2005;85(3):943-978.
3. Foulds IS, Barker AT. Human skin battery potentials and their possible role in wound healing. *Br J Dermatol.* 1983;109(5):515-522.
4. Vanable JW. Integumentary potentials and wound healing. New York: Alan R. Liss;1989. pp:171-224.
5. Zhao M. Electrical fields in wound healing—an overriding signal that directs cell migration. *Semin Cell Dev Biol.* 2009;20(6):674-682. doi: 10.1016/j.semcdb.2008.12.009.
6. Park SS, Kim H, Makin IR, Skiba JB, Izadjoo MJ. Measurement of microelectric potentials for a bioelectrically-active wound care device in the presence of bacteria. *J Wound Care.* 2015;24(1):23-33. doi: 10.12968/jowc.2015.24.1.23.
7. Banerjee J, Ghatak P, Roy S, et al. Improvement of human keratinocyte migration by a redox active bioelectric dressing. *PLoS One.* 2014 Mar 3;9(3):e89239. doi: 10.1371/journal.pone.0089239.
8. Kim H, Makin I, Skiba J, et al. Antibacterial efficacy testing of a bioelectric wound dressing against clinical wound pathogens. *Open Microbiol J.* 2014 21;8:15-21. doi: 10.2174/1874285801408010015.
9. Kim H, Izadjoo MJ. Antibiofilm efficacy evaluation of a bioelectric dressing in mono- and multi-species biofilms. *J Wound Care.* 2015;24 Suppl 2:S10-S14. doi: 10.12968/jowc.2015.24.Sup2.S10.
10. Banerjee J, Ghatak PD, Roy S, et al. Silver-zinc redox-coupled electroceutical wound dressing disrupts bacterial biofilm. *PLoS One.* 2015;10(3):e0119531. doi: 10.1371/journal.pone.0119531. eCollection 2015.
11. 12. Harding AC, Gil J, Valdes J, Solis M, Davis SC. Efficacy of a bio-electric dressing in healing deep, partial-thickness wounds using a porcine model. *Ostomy Wound Manage.* 2012;58(9):50-55.
12. Blount AL, Foster S, Rapp DA, Wilcox R. The use of bioelectric dressings in skin graft harvest sites: a prospective case series. *J Burn Care Res.* 2012;33(3):354-357. doi: 10.1097/BCR.0b013e31823356e4.
13. Whitcomb E, Monroe N, Hope-Higman J, Campbell P. Demonstration of a microcurrent-generating wound care device for wound healing within a rehabilitation center patient population. *J Am Coll Clin Wound. Spec.* 2013;4(2):32-39. doi: 10.1016/j.jccw.2013.07.001. eCollection 2012 Jun.
14. Varhus JD. A novel bioelectric device enhances wound healing: an equine case series. *J Equine Vet Sci.* 2013;34(3):421-430. 10.1016/j.jevs.2013.07.011.

## Ordering Information

---

### JumpStart™ Contact Layer Dressing

Dressing Size (in)	Dressing Size (cm)	Qty/Box	Item Number
1 × 1 Fenestrated	2.5 × 2.5 Fenestrated	10	ABS-4001
1.5 × 8	3.5 × 20	10	ABS-4005
1.5 × 10	3.5 × 25	10	ABS-4006
2 × 2	5 × 5	10	ABS-4002
2 × 5	5 × 12.5	10	ABS-4025
3 × 3	7.5 × 7.5	10	ABS-4003
4 × 4	10 × 10	10	ABS-4004
8 × 8	20 × 20	1	ABS-4008
12 × 12	30.5 × 30.5	1	ABS-4012

Products advertised in this technique guide may not be available in all countries. For information on availability, please contact Arthrex Customer Service or your local Arthrex representative.



This is not veterinary advice and Arthrex recommends that veterinarians be trained in the use of any particular product before using it in surgery. A veterinarian must always rely on his or her own professional clinical judgment when deciding whether to use a particular product. A veterinarian must always refer to the package insert, product label and/or instructions for use before using any Arthrex product. Products may not be available in all markets because product availability is subject to the regulatory and/or veterinary practices in individual markets. Please contact your Arthrex representative if you have questions about availability of products in your area.

Inspired by the body. Powered by electricity. Energized by results.® Patented V.Dox™ technology is owned by Vomaris Wound Care, Inc.

© 2019 Arthrex, Inc. All rights reserved. | [www.arthrexvetsystems.com](http://www.arthrexvetsystems.com) | VLB1-00021-EN\_C